




Peripheral Nerve Bridging Materials: Using Clinical Evidence to Optimize Outcomes



revolutionizing the science of nerve repair™


conduits were designed as a coaptation aid, not as a bridging material

Conduits are hollow tubes, relying on the body to form a fibrin cable.

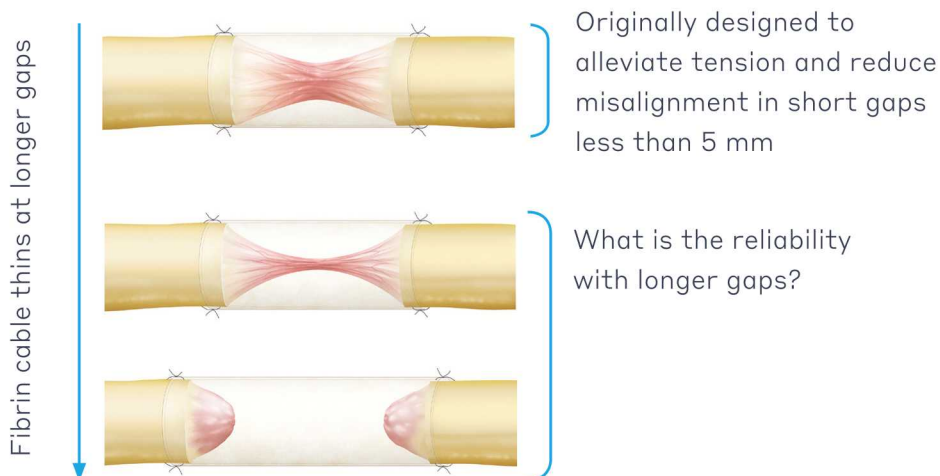

conduits have a reliability threshold

Each datapoint represents one clinical publication in sensory nerve gaps.¹⁻¹⁰

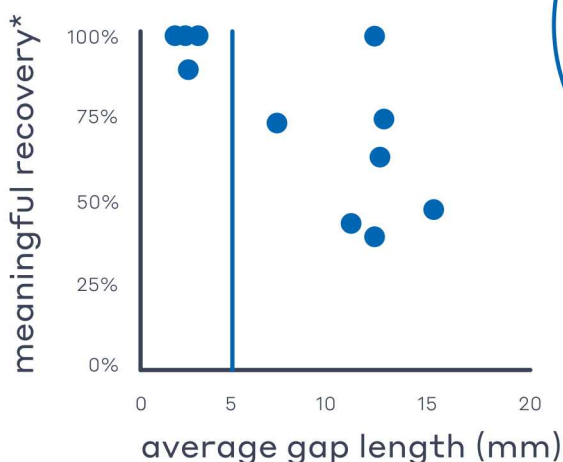

In gaps above 5 mm, more than 1 in 3 conduits fail to achieve meaningful recovery.¹⁻¹⁰

failures may result in:


- Lack of functional recovery
- Absence of protective sensation
- Revision procedure
- Chronic pain

*Defined as S3 or higher.


put conduits back in their place as a coaptation aid in gaps 0–5 mm

avance[®] nerve graft was developed to bridge nerve gaps

It is a biologically active processed human nerve allograft that eliminates the comorbidities and operative time associated with a second surgical site.

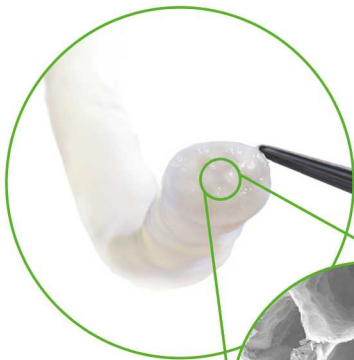


Handles, sutures and flexes at joints similar to native nerve

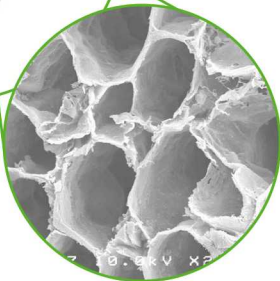

Laminin-rich endoneurial tubes guide regeneration

avance has demonstrated consistency and reliability

Each datapoint represents one clinical publication in sensory nerve gaps,¹¹⁻²⁰

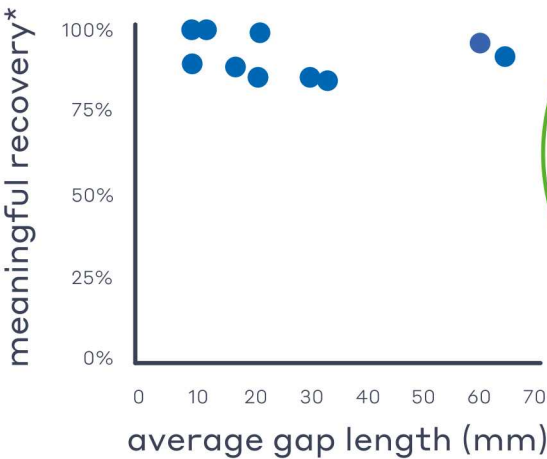

survey says...

Avance is the preferred choice in digital nerve repair over conduit, autograft or vein graft.²¹

results of a survey of 461 hand surgeons with their method for repairing a 2 cm digital nerve gap

*Defined as S3 or higher.

summary of clinical outcomes

Evidence based algorithm avoids potential complications of autograft harvest, without sacrificing patient outcomes.



*Defined as S3 or higher. Data points derived from a weighted average of publications reporting meaningful recovery in sensory nerve gaps.^{1-20,22-24}



references

1. Weber, et al. *Plast Reconstr Surg*. 2000 Oct;106(5).
2. Lohmeyer, et al. *J Reconstr Microsurg*. 2009 Jan;25(1).
3. Alligand-Perrin et al. *Orthop Traumatol Surg Res*. 2011 Jun;97(suppl 4).
4. Taras et al. *J Hand Surg Am*. 2011 Sep;36(9).
5. Chiriac, et al. *J Hand Surg Eur Vol*. 2012 May;37(4).
6. Haug, et al. *J Hand Surg Am*. 2013 Jan;38(1).
7. Lohmeyer, et al. *J Reconstr Microsurg*. 2014 May;30(4).
8. Fujita, et al. *J Oral Maxillofac Surg*. 2014 Jul;72(7).
9. Arnaout, et al. *Chir Main*. 2014 Sep;33(4).
10. Buncke, et al. *J Hand Surg* 2015;40(9).
11. Karabekmez, et al. *Hand (NY)*. 2009 Sep;4(3).
12. Brooks, et al. *Microsurgery*. 2012 Jan;32(1).
13. Cho, et al. *J Hand Surg Am*. 2012 Nov;37(11).
14. Taras, et al. *J Hand Surg Am*. 2013 Oct;38(10).
15. Zuniga. *J Oral Maxillofac Surg*. 2015 Apr;73(4).
16. Rinker, et al. *J Reconstr Microsurg*. 2015 Jun;31(5).
17. Means, et al. *Hand (NY)*. 2016 Jun;11(2).
18. Salomon, et al. *J Oral Maxillofac Surg*. 2016 Dec;74(12).
19. Zuniga, et al. *J Oral Maxillofac Surg*. 2017 Dec;75(12).
20. Rinker, et al. *Ann Plast Surg*. 2017 Jun;78(6S suppl 5).
21. Azouz, et al. *Plast Reconstr Surg Glob Open*. 2018 Aug 6;6(8).
22. Stang, et al. *Arch Orthop Trauma Surg*. 2013 Jun;133(6).
23. Manoli, et al. *Microsurgery*. 2014 Nov;34(8).
24. Kallio. *J Hand Surg. British and European Volume*. 1993;18B.
25. Axogen data on file.

axogen's portfolio of products

options for no transection

1.

A diagram showing a yellow nerve segment with a white, cylindrical Axoguard nerve protector wrapped around it.

axoguard®
nerve protector
2.

A diagram showing a yellow nerve segment with a white, rectangular Avive soft tissue membrane wrapped around it.

avive®
soft tissue
membrane

options for 0 mm to 5 mm

1.

A diagram showing a yellow nerve segment with a white, cylindrical Axoguard nerve connector wrapped around it.

A diagram showing a yellow nerve segment with a white, cylindrical Axoguard nerve connector wrapped around it.

axoguard®
nerve connector
2.

A diagram showing a yellow nerve segment with a white, cylindrical Axoguard nerve connector wrapped around it.

options for 5 mm to 70 mm

1.

A diagram showing a yellow nerve segment with a white, cylindrical Avance nerve graft wrapped around it.

avance®
nerve graft
2.

A diagram showing a yellow nerve segment with a white, cylindrical Avance nerve graft wrapped around it, and a white, cylindrical Axoguard nerve connector wrapped around the graft.

avance®
nerve graft

+

axoguard®
nerve connector
3.

A diagram showing a yellow nerve segment with a white, cylindrical Avance nerve graft wrapped around it, and a white, cylindrical Axoguard nerve protector wrapped around the graft.

avance®
nerve graft

+

axoguard®
nerve protector

options for 70 mm+

1.

A diagram showing a yellow nerve segment with a white, cylindrical Autograft wrapped around it, and a white, cylindrical Axoguard nerve connector wrapped around the graft.

autograft

+

axoguard®
nerve connector
2.

A diagram showing a yellow nerve segment with a white, cylindrical Autograft wrapped around it, and a white, cylindrical Axoguard nerve protector wrapped around the graft.

autograft

+

axoguard®
nerve protector

this is the new normal

At least 66 portfolio peer-reviewed publications

Over 2/3 of US hand and microsurgery fellows were trained on this algorithm in 2018²⁵

Indications and Trademark Disclaimers

Avance Nerve Graft

REGULATORY CLASSIFICATION: Avance Nerve Graft is processed and distributed in accordance with US FDA requirements for Human Cellular and Tissue-based Products (HCT/P) under 21 CFR Part 1271 regulations, US State regulations and the guidelines of the American Association of Tissue Banks (AATB). Additionally, international regulations are followed as appropriate. Avance Nerve Graft is to be dispensed only by or on the order of a licensed physician.

INDICATIONS FOR USE: Avance Nerve Graft is processed nerve allograft (human) intended for the surgical repair of peripheral nerve discontinuities to support regeneration across the defect.

CONTRAINDICATIONS: Avance Nerve Graft is contraindicated for use in any patient in whom soft tissue implants are contraindicated. This includes any pathology that would limit the blood supply and compromise healing or evidence of a current infection.

Axoguard Nerve Connector

INDICATIONS FOR USE: The Axoguard Nerve Connector is indicated for the repair of peripheral nerve discontinuities with gaps up to 5 mm. The Axoguard Nerve Connector is supplied sterile and is intended for single use.

CONTRAINDICATIONS: The Axoguard Nerve Connector is derived from a porcine source and should not be used for patients with known sensitivity to porcine material.

Axoguard Nerve Protector

INDICATIONS FOR USE: The Axoguard Nerve Protector is indicated for the repair of peripheral nerve injuries where there is no gap. The Axoguard Nerve Protector is supplied sterile and is intended for one-time use.

CONTRAINDICATIONS: The Axoguard Nerve Protector is derived from a porcine source and should not be used for patients with known sensitivity to porcine material.

Avive Soft Tissue Membrane

REGULATORY CLASSIFICATION: Avive Soft Tissue Membrane is processed and distributed in accordance with US FDA requirements for Human Cellular and Tissue-based Products (HCT/P) under 21 CFR Part 1271 regulations, US State regulations and the guidelines of the American Association of Tissue Banks (AATB). Additionally, international regulations are followed as appropriate. Avive Soft Tissue Membrane is to be dispensed only by or on the order of a licensed physician.

INDICATIONS FOR USE: Avive Soft Tissue Membrane is processed umbilical cord intended for homologous use as a soft tissue covering.

CONTRAINDICATIONS: Avive Soft Tissue Membrane is contraindicated for use in any patient in whom soft tissue implants are contraindicated.

Disclaimer: Not all products are available internationally.

© 2020 Axogen Corporation. Avance Nerve Graft, Axoguard Nerve Protector, Axoguard Nerve Connector, Avive Soft Tissue Membrane and revolutionizing the science of nerve repair are trademarks of Axogen Corporation. Axoguard Nerve Connector and Axoguard Nerve Protector are manufactured in the United States by Cook Biotech Incorporated, West Lafayette, Indiana.

Axogen Corporation

Phone 888.Axogen1 (888.296.4361)
Fax 386.462.6801
customer@axogeninc.com
www.axogeninc.com



Intended for distribution
in the European Union

revolutionizing the science of nerve repair™

